


Plasti-med[®]
"supporting healthy life"

NEPHRO FLEX[™]


NEPHROSTOMY BALLOON DILATOR

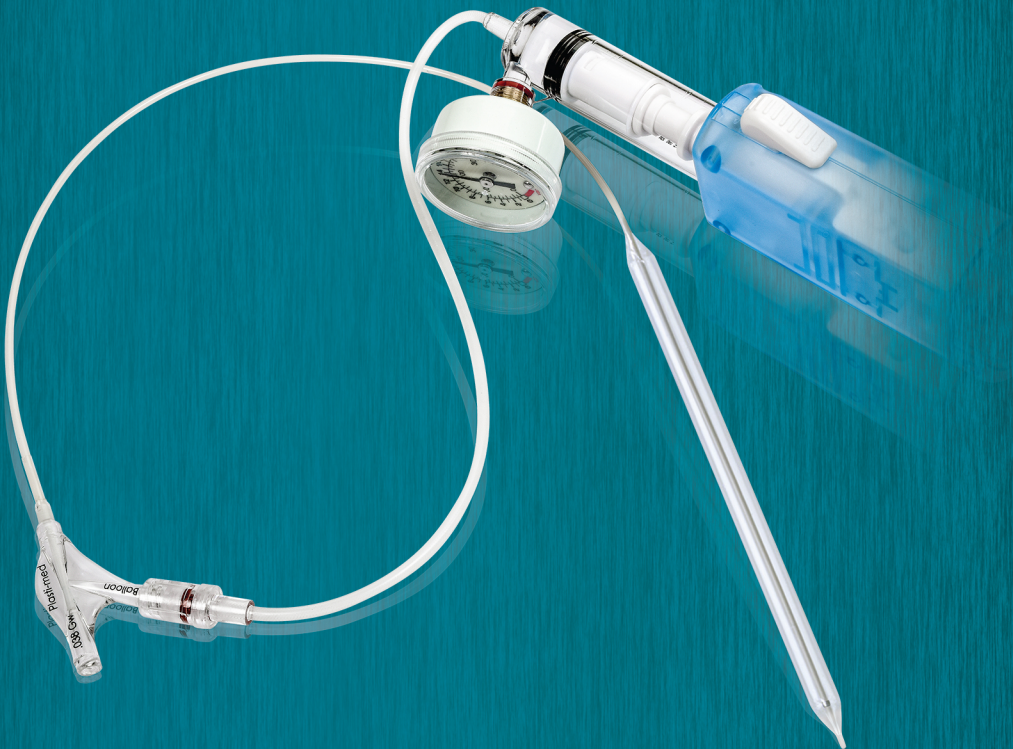

NEPHROFLEX™

NEW ERA IN NEPHROSTOMY WITH LOW TRAUMA AND HIGH EFFICIENCY

NEPHROFLEX™ nephrostomy balloon dilators provides a better approach in endourology surgeries by its safe, fast and less traumatic dilation performance. It comes up with a superior performance even in tightest anatomy and offers accurate size at maximum working pressure.*

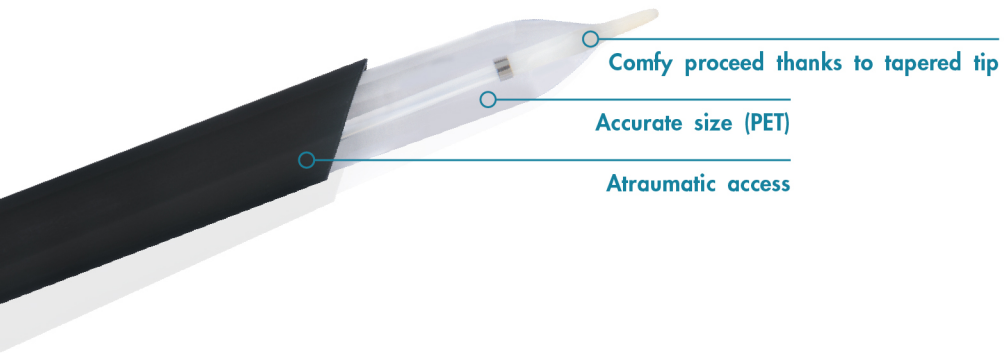

Comfy proceed thanks to tapered tip

Accurate size (PET)

Atraumatic access

- **Platinum/Iridium Marker Band** provides maximum radiopacity to ensure the proper positioning of the balloon in the optimal position within the kidney.

- **Tapered Tip** with angled smooth surface facilitates positioning prior to inflation.


- **Laser and UV welding technology** provides a cleaner and smoother surface than the other bonding methods at the connection points and prevents trauma in the patient.

- **Suitable diameter for 0.038 " guide wire** prevents stuck and strains when reaching the kidney through the wire.

- **Homogeneous inflation and rapid deflation** due to top-notch design of the balloon and high pressure inflation device. This ensures accurate size at every point of the balloon.

*Recommended working pressure is 14 atm and burst pressure is 18 atm.

Inflation Device was designed to provide both rapid inflation and deflation; available in kit and seperately.

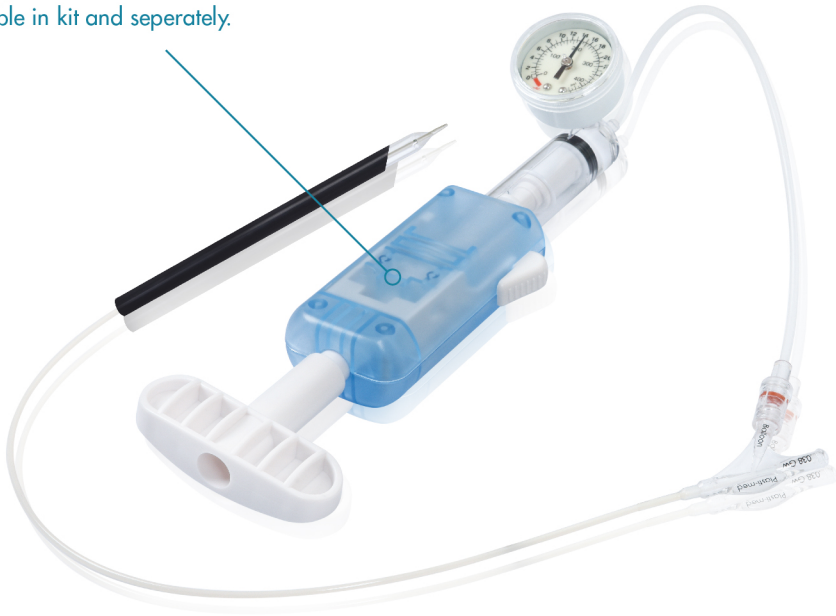

In comparison with traditional method, NEPHROFLEX™ is;

- Safe and faster
- Cost-efficient
- Innovative
- Cause less postoperative complications and provide fast recovery



Code	Name	Balloon Length (cm)	Balloon Inflated (mm/Fr)	Catheter Size (Fr)	Use G.W. (inch)	Usable Length (cm)
533013	Nephrostomy Balloon Dilator	12	10/30	7	0.038	55
533024	Nephrostomy Balloon Dilator with Inflation Device	12	10/30	7	0.038	55
533020	Nephrostomy Balloon Dilator	15	10/30	7	0.038	55
533023	Nephrostomy Balloon Dilator with Inflation Device	15	10/30	7	0.038	55



Since
1995

PLASTİ-MED PLASTİK MEDİKAL ÜRÜNLER SANAYİ VE TİCARET LİMİTED ŞİRKETİ

Deri Osb Mahallesi Yan Sanayi Cad. No: 13 Tuzla/İstanbul, TÜRKİYE

Tel: +90 (216) 591 04 91 **Fax:** +90 (216) 591 04 90

info@plasti-med.com

www.plasti-med.com

